

LINE OF BTI IMPLANTS

Externa®
Interna®



bti.

Human
Technology



CONNECTIONS



The four-lobed **Internal connection** (six-lobed for 3.0 platform) features an absence of protruding edges, which avoids fatigue points. This geometry allows a good load distribution (axial, lateral and torsional) and reduces the risk of plastic deformation during the insertion and de-insertion of the implant. This type of connection also **facilitates the insertion of the prosthodontic components**.

The **external connection** is a hexagon that has different sizes and heights depending on the platform, and are therefore consistent with the biomechanical requirement of each of them.

UNICCA® SURFACE



The **UnicCa® surface** is a chemical modification with calcium ions over its triple roughness

ROUGHNESS 1 · NECK

UnicCa® surface in the coronal area presents an attenuated roughness that enhances marginal tissue retention and reduces bacterial colonization.

ROUGHNESS 2 · THREADS

UnicCa® surface presents a specific porosity on the threads that allows bone anchorage, enhancing primary stability and improving osseointegration.

ROUGHNESS 3 · VALLEYS

UnicCa® surface presents an ideal medium roughness on the valleys that guides bone growth between the threads and at the same time maintains the implant's mechanical properties.



BENEFITS

ELECTROPOSITIVE AND CLEAN

Maintains the superhydrophilic properties.

OSTEOGENIC

Induces the formation of bone tissue.


PERI - IMPLANT BONE STABILITY PROMOTER

Reduces implant failure.

MINIMIZES BACTERIAL ADHESION (in vitro studies)*

[*] Eduardo Anitua, Ricardo Tejero, Miguel Angel Pacha-Olivenza, Maria Coronada Fernandez-Calderon, Maria Delgado-Rastrollo, Mari Mar Zalduendo; Wiley Online Library 2017 Jan. DOI: 10.1002/jbm.b.33860

Apex

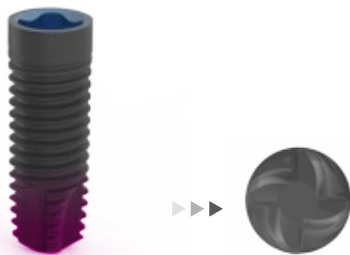

1. Conventional **TAPERED APEX** with great advance capacity, large osseointegration surface and little compression.



2. **CYLINDRICAL APEX** with a lower advance/drive capacity and a large apical surface adapted to achieve primary stability without compressing important anatomical structures (maxillary sinus and dental nerve)

INTERNA®

CONNECTION

3.0 platform

3.0 mm. Ø



Diameter 2.5 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 |

Diameter 3.0 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 |

Diameter 3.3 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 |

Narrow platform

3.5 mm. Ø



Diameter 3.3 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 3.5 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 3.75 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 4.0 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 |

Diameter 4.25 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 |

Diameter 4.75 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 |

Universal platform

4.1 mm. Ø



Diameter 3.3 mm.

Length (mm.) | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 3.5 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 3.75 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 4.0 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 4.25 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Universal Plus platform

4.1 mm. Ø



Diameter 4.5 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 5 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 5.5 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 |

Wide platform

5.5 mm. Ø



Diameter 5.5 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 6 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 |

Diameter 6.25 mm.

Length (mm.) | 5.5 | 6.5 | 7.5 | 8.5 |



EXTERNA[®]

CONNECTION

Tiny[®] implant
Narrow platform

3.5 mm. Ø



Diameter 2.5 mm.

Length | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 3 mm.

Length (mm.) | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 3.3 mm.

Length (mm.) | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 3.5 mm.

Length (mm.) | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 3.75 mm.

Length (mm.) | 4.5 | 5.5 | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Universal platform

4.1 mm. Ø



Diameter 3.75 mm.

Length (mm.) | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 4 mm.

Length (mm.) | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Universal Plus platform

4.1 mm. Ø



Diameter 4.5 mm.

Length (mm.) | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

Diameter 5 mm.

Length (mm.) | 6.5 | 7.5 | 8.5 | 10 | 11.5 | 13 | 15 |

PREDICTABLE IMPLANTOLOGY

THE WIDEST AND MOST VERSATILE RANGE ON THE MARKET

Various prosthodontic platforms depending on the prosthetic solution we want to obtain and various implant sizes depending on the bone substrate we have. **The implant system adapts to the residual bone volume of your patient.**

A GUARANTEE FOR EXTREME RESORPTIONS

TRANSVERSAL ATROPHY

Narrow implants for provoking ridge expansions and increasing the volume of the bone substrate without having to perform less predictable surgeries.

VERTICAL ATROPHY

Short and extra-short implants allowing to perform surgeries in maxillary regions with a lack of alveolar bone height

THE PERFECT TEAM: ENDORET® (PRGF®) + UNICCA® IMPLANTS

The combination of ENDORET technology, the pioneer in obtaining plasma rich in growth factors that stimulates and accelerates tissue healing and regeneration, and the nano-roughness surface of the UnicCa® range of implants **achieves higher osseointegration rates in less time and reduces bacterial colonization.**



CleanImplant Trusted
Quality Mark

2017-2019

UnicCa® surface, certified pureness

BTI UnicCa® is the first implant system in the market awarded with the CleanImplant Foundation Mark, as a guarantee of the highest quality of its materials and surface.



BTI Comercial

San Antonio, 15 - 5º
01005 Vitoria-Gasteiz
(Álava) · SPAIN
Tel: +34 945 140 024
Fax: +34 945 135 203
pedidos@bticomercial.com

B.T.I.

Biotechnology Institute S.L.
Parque Tecnológico de Álava
Leonardo da Vinci ,14
01510 Miñano (Álava) SPAIN
bti.implantes@bti-implant.es

Subsidiaries

GERMANY

Mannheimer Str. 17
75179 Pforzheim · Germany
Tel. +49 7231 428060
Fax +49 7231 4280615
info@bti-implant.de

PORTUGAL

Praça Mouzinho de Albuquerque 113, 5º
4100-359 Porto · Portugal
Tel: (351) 22 120 1373
Fax: (351) 22 120 1311
bti.portugal@bticomercial.com

FRANCE

6 Avenue Neil Armstrong
Immeuble Le Lindbergh
33692 Merignac CEDEX · France
Tel: (33) 556 18 11 18
info@bti-implant.fr

UK

870 The Crescent
Colchester Business Park · Colchester
Essex CO49YQ · UK
Tel: (44) 01206580160
Fax: (44) 01206580161
info@bti-implant.co.uk

ITALY

Via Conservatorio 22
20122 Milano · Italy
Tel: (39) 02 7060 5067
Fax: (39) 02 7063 9876
bti.italia@bti-implant.it

USA

1730 Walton Road
Suite 110
Blue Bell, PA 19422-1802 · USA
Tel: (1) 215 646 4067
Fax: (1) 215 646 4066
info@bti-implant.us

MEXICO

Ejercito Nacional Mexicano 351, 3A
Col. Granada Delegación Miguel Hidalgo
Messico DF · CP 11520 · Mexico
Tel: (52) 55 52502964
Fax: (52) 55 55319327
bti.mexico@bti-implant.com

www.bti-biotechnologyinstitute.com

NOTE: consult your distributor for the availability of the product in the different markets.



BTI APP
NEW PRODUCTS

iPhone / smartphone version
iPad / Tablets version (**Customer Area**)